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Fairywill

P11
Ring the future

Thank you for choosing Fairywill P11 sonic toothbrush. Please read the instructions carefully before using the product, and keep the user manual in a safe place for future reference.

IMPORTANT SAFEGUARDS

1. This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
2. Children should be supervised to ensure that they do not play with the appliance.
3. This appliance contains batteries that are non-replaceable. When the battery is at end of life, the appliance shall be properly disposed of.
4. Do not expose the appliance or battery to excessive temperatures.
5. Be aware of the risk of terminals of the battery-operated appliance or battery being short circuited by metal objects.
6. The appliance is must only be supplied at safety extra low voltage corresponding to the marking on the appliance.
7. The power supply with output power not exceeding 15W in normal and abnormal working conditions must be used.
8. Periodically check the cord for damage. If cord is damaged, take the charging unit to an Service Center. A damaged or non-functioning unit should no longer be used.
9. If the toothbrush handle is dropped, the brush head should be replaced before the next use even if no damage is visible
10. Do not touch the charging line with wet hands. This can cause an electric shock.
11. Do not modify or repair the product. This may cause fire, electric shock or injury.
12. Please use an adapter so that the voltage is 5V==1000mA during charging, otherwise the battery could be damaged by overcharging.
13. If the battery is damaged or leaking, avoid contact with the skin or eyes. If this occurs, immediately rinse well with water and seek medical care.
14. This appliance is just for household use.

SAVE THESE INSTRUCTIONS

There is a built-in rechargeable lithium battery in the product. When disposing of the appliance, do not mix it with general household waste. It should be sent to the battery recycling center for disposal in accordance with local disposal management regulations.

The symbol on the product or on its packaging indicates that this product may not be treated as household waste.

Instead it should be handed over to the applicable collection point for the recycling of electrical and electronic equipment. By ensuring this product is disposed


of correctly, you will help prevent potential negative consequences for the environment and human health, which could otherwise be caused by inappropriate waste handling of this product. For more detailed information about the recycling of this product, please contact your local council office or your household waste disposal service.

1. Features:

- Waterproof and non-slip design.
- Replaceable brush heads.
- High-frequency vibration delivers 62,000 brush movements each cycle*, removing plaque effectively, preventing gingivitis, and cleaning your teeth and gums deeply.
- USB charging, travel friendly and compatible with global voltage.
- A low-power reminding function, when the voltage of the battery is low on power, the indicator light will flash slowly or the toothbrush is unable to turn on during the operating process, as a reminder to charge the toothbrush.
- A 2 minute smart timer with interval pausing, that gives a reminder to change the brushing location every 30 seconds.
- There are 3 modes: Clean, Soft and White settings built into the toothbrush. There is also a mode memory function built into the brush, that means when you turn on the toothbrush, the default starting mode stays the same as the mode that was set at the time the toothbrush was last shutoff.

2. Designs: (Please see pic. 1)

- ① Brush Head
- ② Power Button Mode Button
- ③ Function(Mode)/Power/Charging
- ④ Indicator Light
- ⑤ Charging Port

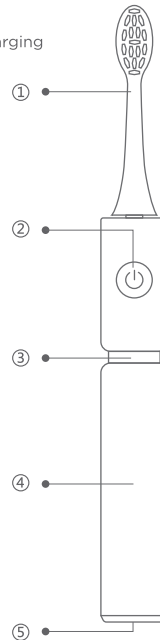pic. 1

3. Operations:

★ Using the toothbrush

- Slide the toothbrush head down onto the handle and make sure it is locked in the correct position.
- Turn on the toothbrush by pressing the “Power” button.
- Shift between the 3 modes (“Clean, Soft and White”) by continuously pressing the “Power” button.
- Wet the bristles and apply an appropriate amount of toothpaste onto the brush.
- Place the bristles with toothpaste against the teeth at a slight angle towards the gum line.
- Gently move the brush head across the teeth on the front, back, and all chewing surfaces. The brush should be allowed to glide across the teeth with minimal pressure. Do not apply too much pressure to the teeth/gums, or apply the bristles to one area of the mouth for an extended period of time.
- The brush has a 2 minute smart timer as this is recommended by dentists as the optimal time for brushing the teeth and there is interval pausing after every 30 seconds as a reminder to move to another quadrant of the mouth.
- Move to the next section when you feel the pause and the toothbrush will automatically power off after 2 minutes.
- To avoid splashing toothpaste, only turn on the toothbrush after you have placed it against your teeth and only remove it from the mouth after it is turned off.

★ Changing the toothbrush head

- Install: Slide toothbrush head down into the handle. Please make sure it is in the correct position.
- Detach: Pull the toothbrush head out along the handle.

★ Charging the toothbrush

- If the toothbrush is in low power when turning it on for the first time, please ensure it is fully charged before usage.
- If the toothbrush is not in use for more than 3 months, please fully charge it according to the charging instructions before using it.
- If the indicator light flashes slowly during the operation process, it indicates the toothbrush is in a low battery condition, please charge the toothbrush as soon as possible.
- When you are going to charge your toothbrush, please insert the USB cable into the charging port at the bottom of the toothbrush and ensure the brush is connected to the electric power source safely.
- The charging indicator light will be white and flashes slowly during the charging process. When the battery is fully charged, the indicator will come on for two minutes and not flash anymore, after 2 minutes, it will turn off automatically. After the toothbrush is fully charged, it is possible to use it for two minutes per section, twice a day for 30 days .
- There is a waterproof function built into the toothbrush, so it can be used in the bathroom safely.

★ Modes Introduction

- The first mode: Clean
In this mode the brush delivers about 62,000 brush movements each cycle*. 3 white indicator lights turn on, the swing of the brush is large and the vibration intensity is the strongest. It is suitable for normal oral conditions, to gain the strongest cleaning effect, and maintain healthy oral hygiene.
- The second mode: Soft
In this mode, the brush delivers about 56,000 brush movements each cycle*, there is only 1 white indicator light that turns on, the brush swing is soft and the vibration intensity is the softest of the three modes, which is suitable for users with sensitive teeth.

- The third mode: White
In this mode, the brush delivers about 56,000 to 62,000 brush movements each cycle*, there are two white indicator lights turned on, the bristle swing and the vibration intensity changes rapidly. For gentle stimulation, designed to keep gums firm and healthy.

4. Specifications:

Working temperature: 0-50 C
Working voltage: DC3.7V
Bristles material: Nylon
Input voltage: DC5.0V
Rated power: 0.7-1.8W
Waterproof rating: IPX7

5. Precautions:

- Do not disassemble, repair, or alter the product without authorization.
- People who have severe periodontal disease or oral surgery should consult a doctor before use.
- Toddlers, infants and people who are unable to operate should be supervised by a parent or adult.
- Do not soak the product in water for an extended period of time.
- Do not place the toothbrush in high temperature condition or expose it to the sun.
- Do not force the bristles hard against the teeth or gums.
- Not suitable for children under 12 years old.

6. After-sales service:

- A 1-year warranty starts on the day of purchase, please contact Fairywill when it is necessary.
- Within the warranty period, any defect arising from materials or product process, we will replace the parts or the whole product depending on the condition.
- When you apply for a free repair, you must submit the toothbrush and proof of the date of purchase (purchasing invoice) to Fairywill.
- The following conditions are not included in the warranty:
 - Damages caused by improper use, for example, improper use of the power supply, plugged into improper socket, artificial damage, etc.
 - Ageing and wearing-out caused by normal use.
 - Damages due to product being used for other purposes, or use of components not supplied by Fairywill.

*A 2 minute brushing cycle is recommended by dentists.

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